

INSTRUMENT STAND Model 2000-ST

RIGHTmed
S-OPTIK



S-OPTIK

FEATURES:

◆ logical INSTRUMENT CONSOLE

durable membrane switches, logically positioned, control power mains, chair up/down and auto-return, overhead lamp with rheostat, and binding post (BIO) control with rheostat. Three rechargeable wells accommodate all brands of hand held instruments.

◆ intelligent SLIT LAMP ARM

an engineering achievement, our fully counterbalanced arm, incorporates the full function and range, with the lock/release mechanism conveniently located at your finger tips. A 13" (33cm) range of travel and counter balance range up to 60lbs (27.2kg), this arm easily accommodates all, including wheel chair patients. The chair up/down switch is conveniently located on the front of the arm and is reversible.

◆ ergonomic UPPER (REFRACTOR) ARM

engineered with natural, and comfortable hand movements in mind, the lock/release lever engages a single lock mechanism in all three dimensions (forward/backward/vertical). A 12.5" (31.75cm) vertical range, and counter balance up to 20lbs (9.1kg) ensures precise and effortless positioning.

◆ elegant HALOGEN OVERHEAD LIGHT

the overhead light completes the modern appearance of the 2000-ST, yet offers a long-lasting halogen lamp, providing bright, direct illumination for general lighting or focused minor procedures. Metal lamp arm construction allows stable, one handed movements with superior durability.

◆ optional THIRD ARM


with the same ergonomically designed lever mechanism as the refractor arm, the optional third arm delivers a vertical range of 12.5" (31.75cm) and counter balance up to 30lbs (13.6kg) to easily accommodate the addition of a Keratometer to your examination room.

◆ sensible OPTIONAL ADDITIONS

STORAGE TRAY



WHEELCHAIR KIT

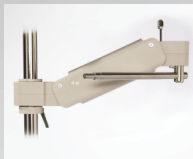

BIO HANGER



THIRD ARM



SPECIFICATIONS

Measurements

| | |
|-----------------|------------------------------------------------------------------|
| Base | 22" (55.8 cm) wide x 23.25" (59.6 cm) depth |
| Height | 34.5" (87.6 cm) to top of console, 68" (172.7 cm) to top of pole |
| Shipping Weight | 240 lbs (108.8 kg) gross 255 lbs (102.1 kg) net. |

TECHNOLOGY
VALUE
 RELIABILITY

Distributed by:

with a policy of continuous improvement, we reserve the right to make changes to specifications without advance notice.



5325 Cleveland Street, Suite 303 Virginia Beach, VA 23462
 1.888.224.6012 | www.s4optik.com