

5300

Stool Series

800.735.0357
haag-streit-usa.com


Model 5340
Seat Substrate 14-Ply Birch



Exceptional Stability from Five-legged Design

Backrest Adjusts Horizontally and Vertically (Models 5346 And 5348)

Pneumatic Lift

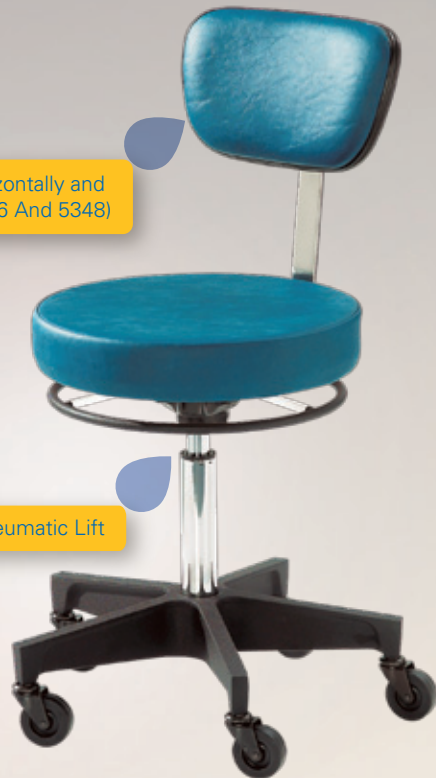

Ergonomic Wrap-around Arm

Model 5348
Comfort and Support From High-density Polyfoam Waterfall Seat



Model 5356
Easily Adjust Height with Control Ring



53000

Stool Series



800.735.0357
haag-streit-usa.com

◆ **Exceptional Stability**

◆ **Finish and Upholstery**

- Select from 14, easy-to-maintain vinyl colors. Other vinyl colors available by request.
- All exposed metal with exception of the aluminum base is chrome-plated steel.
- Rugged 2½" hard rubber casters.

◆ **Easily Adjustable Height Mechanism**

Lift the control ring located under the seat. The stool's pneumatic device will help lift the seat as you ease your weight off the stool. Vertical travel is 6.25". Seat height ranges from 18.50" to 24.75". Footrest ring available as an option.

◆ **Tech Specs**

Base diameter 19¾". Removable footrest: 16" diameter. Shipping weights: 5340 (26 lbs.), 5346 (32 lbs.), 5348 (32 lbs.), 5356 (47 lbs.)

	5340	5346	5348	5356
Seat Dimensions	15"	15"	16 x 14" sloped waterfall construction supports thighs	15"
Seat Construction	2" polyfoam	2" polyfoam	2" polyfoam	2" polyfoam
Swivel	360°	360°	360°	360°
Five legs for extra stability	Yes	Yes	Yes	Yes
Circular control ring for pneumatic lift	Yes	Yes	Yes	Yes
Upholstered 3-way, wrap around arm	No	No	No	Yes, exceptional abdominal, back or arm support

Examination Stool with Pneumatic Assisted Lift.